270mm

135mm

\_teknofit

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EMS 4010

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
EMS 4010

Calf Relaxation

Thank you for choosing our massager.

Please read the user manual carefully before use.
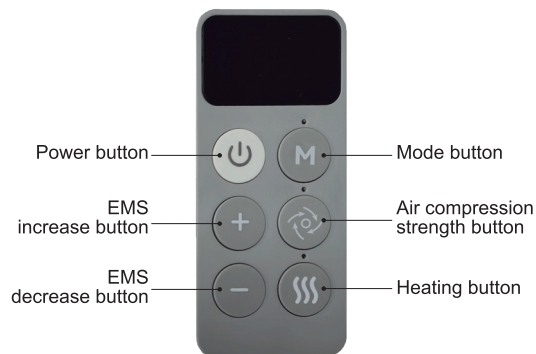
#### HOW TO OPERATE THE CONTROLLER

The device has 5 modes, 9 adjusted levels of EMS intensity, 3 strengths of air compression, heating function and works for 15 minutes.

Press the power button to switch on the device. The screen displays "--".

Press the power button to switch off the device. The screen will shut down.

After switching on the device, press the "mode" button to select mode, totally five modes.

After selecting mode, press the "+" or "-" button to adjust the EMS intensity. There are 9 levels of output intensity. When the intensity is adjusted, the screen will display level of intensity accordingly.

Press the air compression strength button to adjust strengths, totally three strengths. When the air compression strength is adjusted, the screen will display levels of strength accordingly.

When the device is on, press the "heating" button to turn on/ turn off the heating function. There is only one level of heating (38°C).

#### CHARGING METHOD

Please use the attached USB cable to charge the device before use. The USB cable can be connected to a computer or adapter and ensures that each charge is enough for 2.5 hours.

Screen display status during charging: the battery signal will flash when the device is charging and will stay on once fully charged.

The device can be used continuously for about 2 hours after being fully charged.